

Trigger Adapter for Ocean Optics Spectrometers

Access the trigger signal from your spectrometer through a BNC connector for coaxial cable

Soluciones y Tecnologías de Control Embebido S.A.P.I. de C.V.




www.sepradel.com

contacto@sepradel.com

cel +52 1 833-3894054

fijo +52 833 116 3738


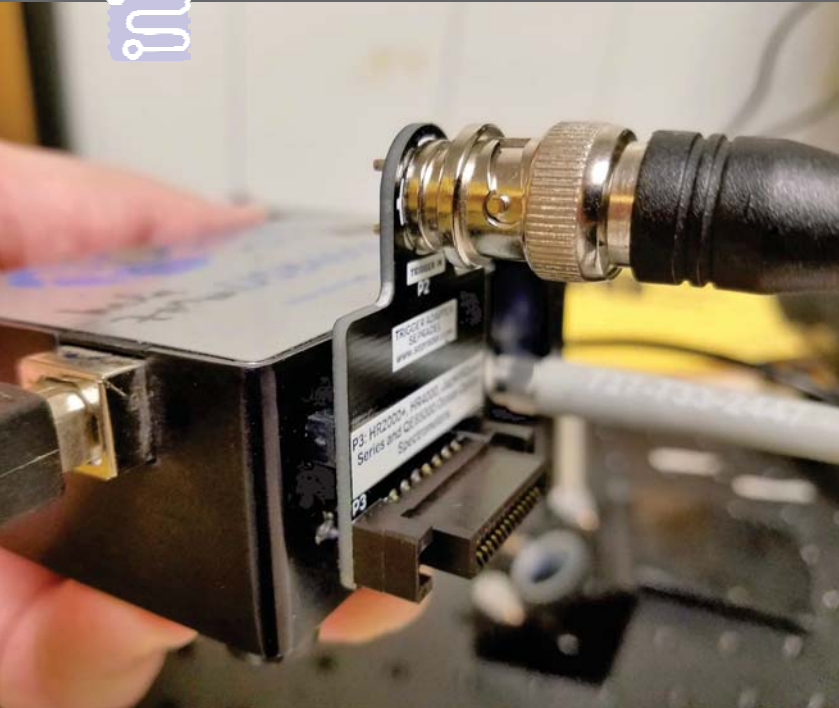


Our adapter, model OOTRIGGER, brand Sepradel® allows access to Trigger Signal of Ocean Optics Spectrometers; While providing a BNC connector for a standard coaxial line. The Sepradel® Trigger Adapter is compatible with the following Ocean Optics spectrometers:

USB2000+,
USB4000,
Jaz
Hr2000

Hr2000+,
Hr4000,
NIRQuest Series,
Qe65000

Connecting the Adapter according to your Model of Spectrometer

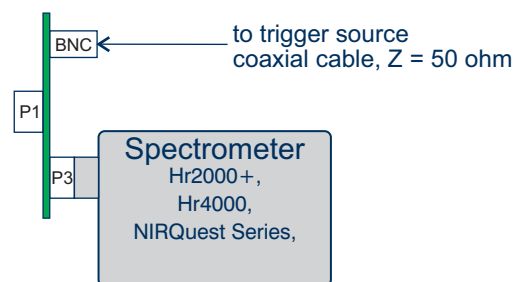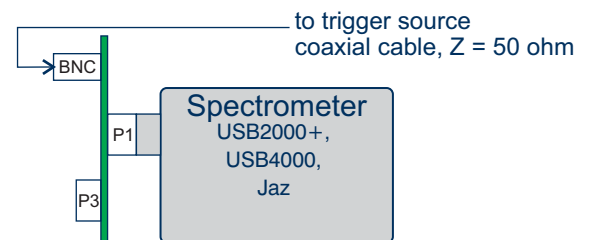

Technical characteristics

- Connection to the external trigger source with BNC connector.
- It Includes RC lowpass filter, R = 50 ohm and C = 330pF.
- Maximum dimensions of 32 mm x 45 mm (width x height).
- Supported models: USB2000+, USB4000, Jaz, HR2000, HR2000+, HR4000, NIRQuest Series and QE65000.

Nota: Ocean Optics, USB2000+, USB4000, Jaz, HR2000, Hr2000+, HR4000, NIRQuest Series and QE65000 Spectrometers are trademarks of Ocean Optics. Neither Sepradel™, nor any software programs or other goods or services offered by Sepradel™, are affiliated with, endorsed by, or sponsored by Ocean Optics.